


OUR TRUSTED LEGACY IS YOUR SUCCESS



PRODUCT CATEGORY:
HOSPITAL FURNITURE

PRODUCT NAME:
HU-I-U0002153

PRODUCT DESCRIPTION:
MOBILE INFUSION BAG HOLDER

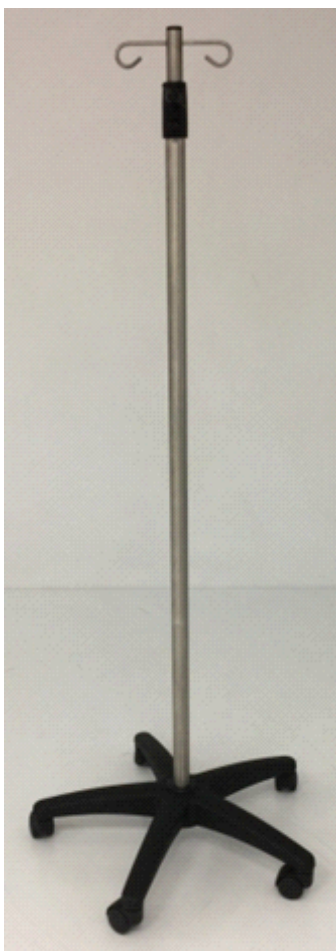

Area of Use: All Areas

Intended Use: The mobile infusion stand is intended to support and suspend intravenous fluid bags, blood bags, or other infusion containers during intravenous therapy and blood transfusion procedures. It is designed for use in healthcare facilities to facilitate the safe and efficient administration of fluids to patients. The stand is mobile, allowing easy positioning close to the patient.

OUR TRUSTED LEGACY IS YOUR SUCCESS



HU-I-U0002153 SPECIFICATIONS

TECHNICAL SPECIFICATION:

Movable height adjustable infusion stand.

Heavy carriage mounted on 5 swivel castors.

Lower - end support column is deep and securely fixed into the carriage - base.

Solid manual lever at the upper - end of the support column, allows setting telescopic upper part at the required height.

A brake prevents exceeding the maximum height setting.

Double hook fixed at the top of telescopic rod.

MATERIALS:

High resistance to corrosion (tropical environment).

Support and telescopic column:	Austenitic stainless steel 18/10 (304 Quality).
--------------------------------	---

Carriage - base:	Single molded unit, polypropylene.
------------------	------------------------------------

Caster frame / bracket:	Polypropylene.
-------------------------	----------------

Wheel bearing:	Polypropylene.
----------------	----------------

Castor fixation into carriage base:	Stainless steel.
-------------------------------------	------------------

All furniture components contain less than 100 parts per million of the below chemical groups, as applicable:

Urea formaldehyde.

Heavy metals including mercury, cadmium, lead, antimony.

Hexavalent chromium in plated finishes.

Stain and non-stick treatments derived from Perfluorinated Compounds (PFCs), including Perfluorooctanoic Acid (PFOA).

Added antimicrobial treatments.

DIMENSIONS:

Height, adjustable:	135 - 225 cm.
---------------------	---------------

Carriage-base, diameter:	60 cm.
--------------------------	--------

Support column:	3.2 cm (outside, across), 1.5 mm (thickness)
-----------------	--

Telescopic upper part:	2.5 cm (outside, across), 1.5 mm (thickness)
------------------------	--

Swivel castor twin - wheel:	4.2 x 5.0 cm (w x diameter) (Plus and minus %5).
-----------------------------	--

Carrying capacity:	10 kg
--------------------	-------

Knockdown construction:	Yes.
-------------------------	------

OUR TRUSTED LEGACY IS YOUR SUCCESS



HU-I-U0002153 SPECIFICATIONS

SUPPLIED WITH:

- 1 x complete set of tools required for assembly.
- List of accessories and parts.
- Detailed step-by-step instructions for assembly and safe use, text-free pictorial based (i.e. line-drawings only).
- Instructions for use including cleaning and disinfection instructions in English, French and Spanish.

PACKAGING, LABELLING, INSTRUCTIONS:

- One (1) unit per box
- Unit Weight in 7.1 Kg (including its packaging)
- Unit Volume in 0.12 m3 (including its packaging)
- Symbols in accordance with ISO 15223
- The packaging made of recycled or recyclable materials.