



PRODUCT CATEGORY:
HOSPITAL FURNITURE

PRODUCT NAME:
HU-R-011

PRODUCT DESCRIPTION:
STANDARD DIALYSIS & CHEMOTHERAPY CHAIR



Area of Use: Renal Units / Cancer Wards / Examination Rooms

Intended Use: Exam Chair is designed to provide positioning and comfortable support for patients during general chemotherapy and dialysis examinations and procedures by qualified caregivers.



SAFETY FEATURES

ARMREST MOVEMENT FLEXIBILITY

Armrest can be adjusted to move up and down by using plunger and getting on or off the recliner easily.



WHEN LYING ON THE RECLINER:

The armrest is on standby position to act like side-rails for patient safety, to reduce the risk of falling.



WIDE ARMREST AREA:

Allowing more treatment's area to the nurse so both patient and nurse can feel more comfortable during treatment process.

EMERGENCY BACK LOWERING LEVER

This function is to lower the backrest up to -5° so it can be used rapidly in case of emergency, such as when patient is fainting due to anaemia or similar conditions or even a power failure.




PULLING THE LEVER :

When pulling the lever, hold the bed to ensure smooth movement while adjusting the seat for the patient.

LEVER POSITION:

Available on one side of the head rest edge.

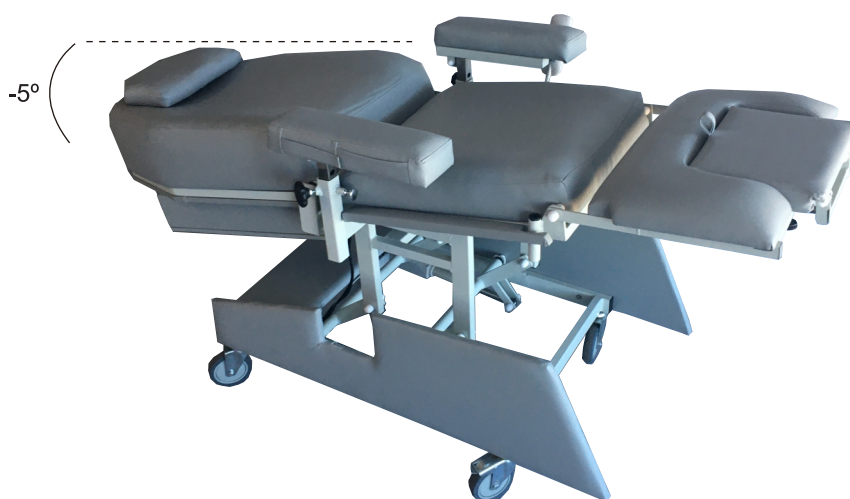

EMERGENCY BACK LOWERING LEVER MAXIMUM ANGLE:

Maximum back raise angle by using emergency lever.



COMFORT FOR PATIENTS

WIDER SEAT SIZE :

The wide space seat is designed to enhance the comfort of the patient while taking treatments and to have more range of movement.



ACCESSORIES POCKET:

Designed for simple storage of patient's belonging. Example : handphone, keys, pen, and other things. This pocket have an opening at its bottom for an easier cleaning process and cable lane while charging your handphone or similar needs with USB port.



PULLING THE LEVER :

There are two holes on the frame side that can be used to adjust footrest position according to patient height by using knob bolt.



ADJUSTABLE PILLOW HEIGHT:

To make it easier for the patient to adjust the position of the pillow with one hand, a free-stop pillow with a weight for balancing is used, and the pillow is attached by a belt carrying the weight.






WORK EFFICIENCY

HYGIENE SEAT COVER SURFACE :

In the interest of hygiene, the seat cover is made with a seamless material that is impermeable. To clean soiled seat covers, simply wipe away the contaminating liquid.

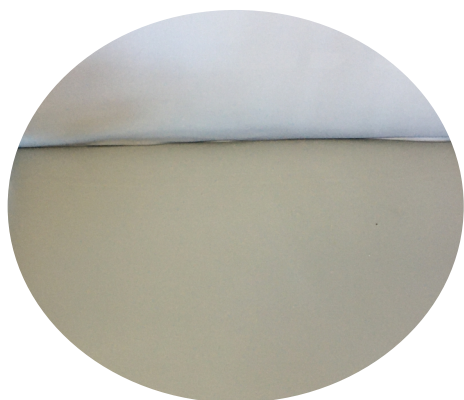

HANDHELD CONTROLLER :

Handheld controller features a slim width ergonomically sculpted grip, making it easier to clean.



ADDITIONAL ACCESSORIES



IV POLE



PILLOW



HU-R-011 SPECIFICATION

SPECIFICATIONS:	
Motors:	3 Motors
Total Length:	1270 mm (When flat 2140 mm and with the foot piece extended)
Total Width:	830 mm
Total Height (*1):	1340 mm
Seat Width:	590 mm
Seat Height:	570 mm
Castors:	125 mm diameter castors, 2 with brakes.
Power Cable Length:	2000 mm
User's Maximum Weight:	160 kg
Safe Working Load:	210 kg
Back Raise Angle (*2):	0° to 75°
Leg Raise Angle:	-30° to 0°
Mattress Base Material:	Particle Board with Melamine
Main Frame Material:	Made of Steel, Powder Coating
Arm Rest Material:	Vinyl Leather, Urethane Foam



TECHNICAL DRAWINGS:

